

Kiosk Pro Lite Set Up



To collect patient data on iPads in clinic securely Amplitude recommends using a combination of Kiosk Pro Lite and Guided Access. Kiosk Pro Lite can be installed for free in the Apple App Store (IOS 18 or later)

1. Open the Apple App Store on the iPad.

2. Search for and download the Kiosk Pro Lite app.

3. While the app is downloading open the iPad's Settings.





4. Select General and then Accessibility.

5. Scroll down and select Guided Access.

6. Switch Guided Access on and set a passcode to switch this off when needed.

7. Close the iPad Settings.

8. To set up Kiosk Pro Lite for Amplitude you will need to obtain the link to the Patient Portal. (See next page)

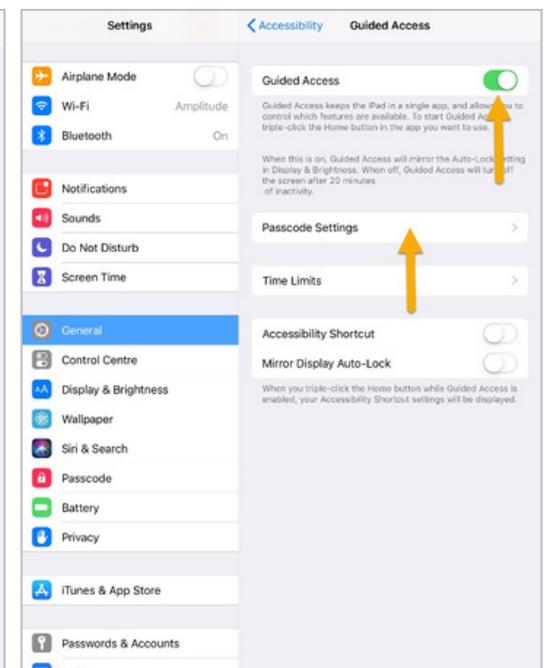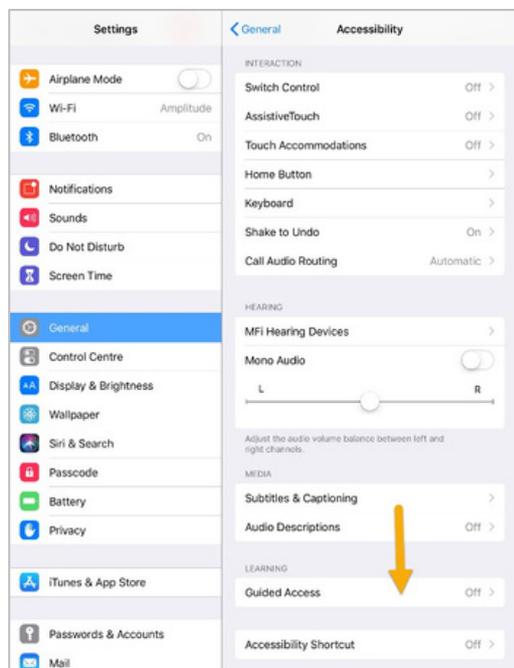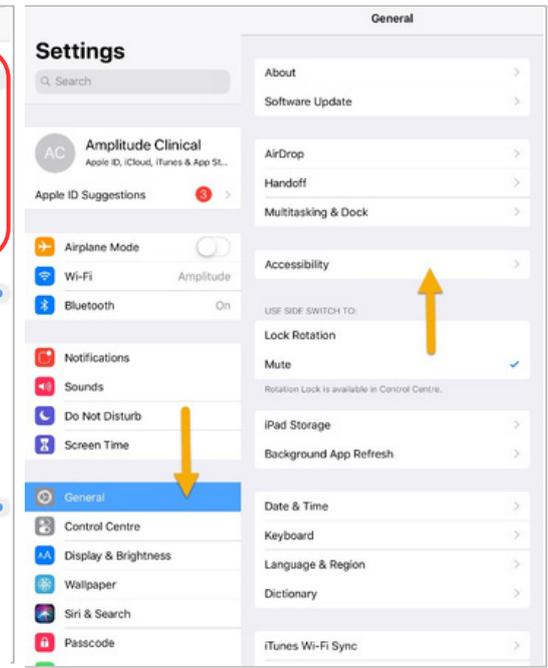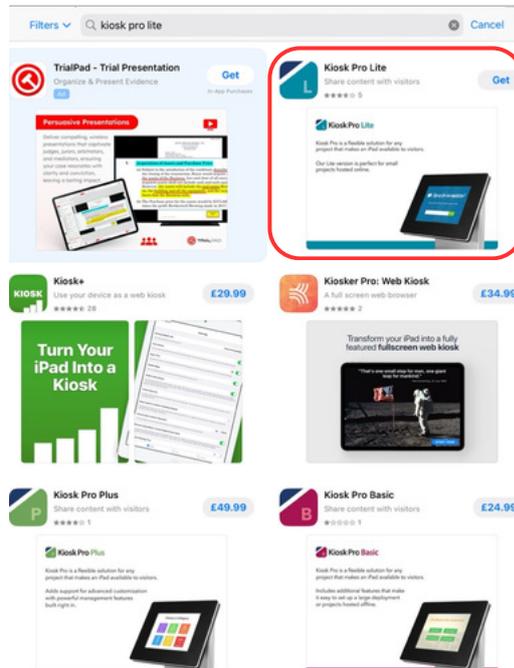

9. To obtain this open Amplitude and click Settings. Click System Settings and copy the Patient Portal URL.

10. Paste the URL into the Homepage Content'.

11. Launch the Kiosk Pro Lite app.

12. Start Guided Access by triple-clicking the physical Home button (for iPads without a home button, triple click the Power button).



13. Under 'Options' you can choose to disable 'Sleep/Wake Button' if you want to prevent the device from being put to sleep using the power button.

14. Tap 'Start'.

12. Run the Kiosk Presentation

13. This should display the Patient Portal page.

14. Triple click the iPad Home Button to activate Guided Access, this stops any patients using other apps on the iPads in clinic.

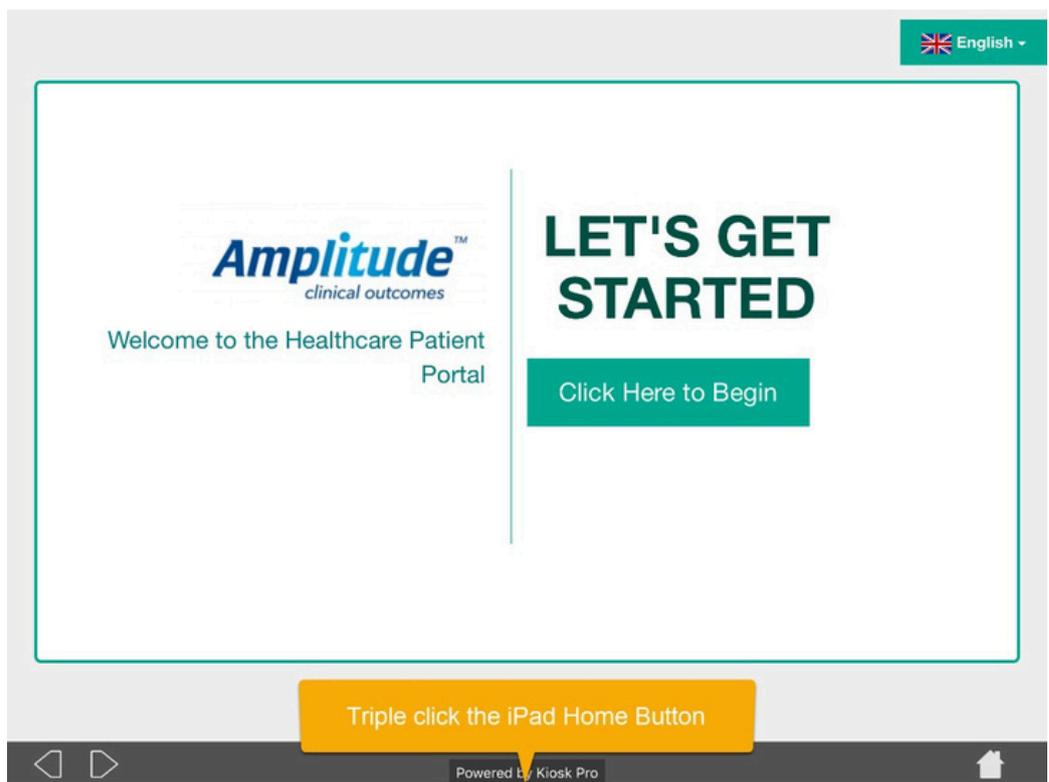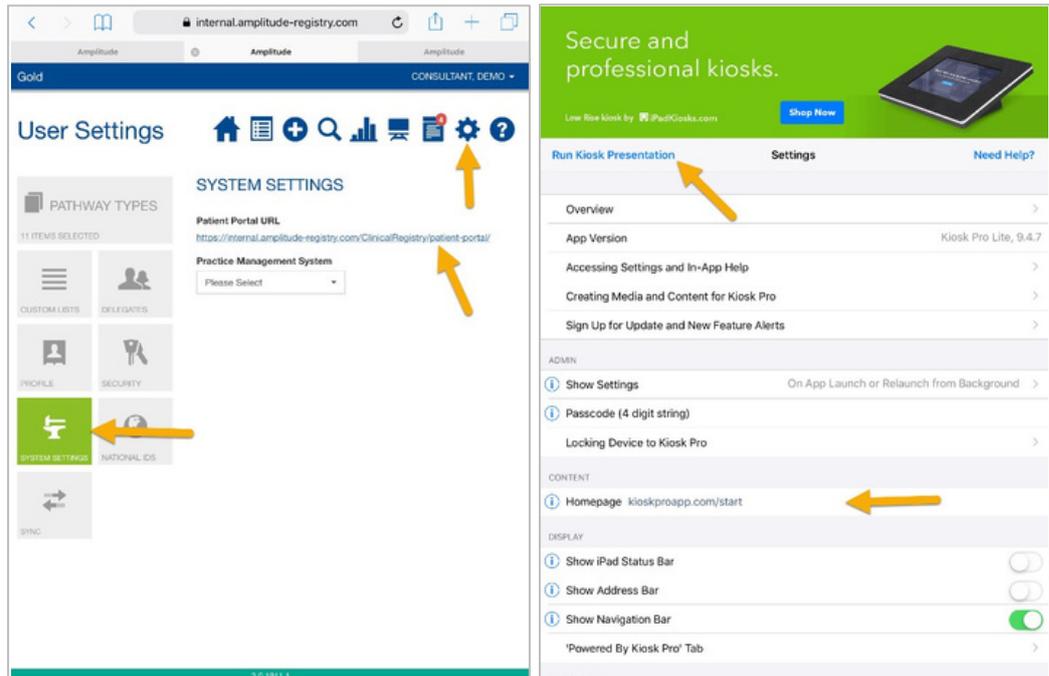

To exit Guided Access, triple-click the Home button (or Power button) and enter the passcode.

The on-screen notifications when a user touches the Home button cannot be turned off and may not work for your intended kiosk presentation. As this is part of the operating system, it is not something we can change from within the app.